

NUMBERED NOTICE


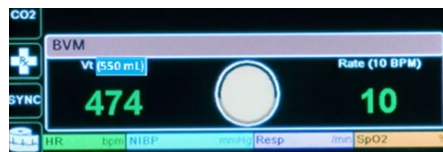
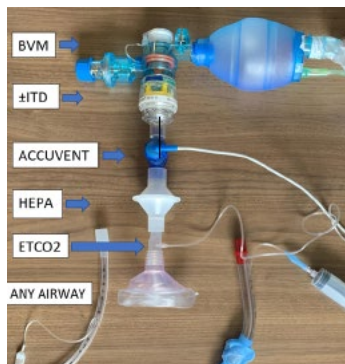
Number: **23-01**
Date: 4/24/2023
From: Ben Weston, MD, MPH, Medical Director
To: All EMS Provider
Subject: Accuvent Clinical Trial
Purpose: Directive

Thank you to all who have participated so far in this clinical research trial to better understand outcomes in adult out of hospital cardiac arrest. The Zoll Monitor Real BVM HELP Phase II begins over the week of April 25th.

Reminders: **DO NOT USE ON CHILDREN**

Once the BVM ACCUVENT Help feedback screen is turned on, remember:

1. Apply to adults only (presumed or known to be ≥ 18 years old)
2. Target ventilation volume goal is 550 ml per breath and the device is programmed for this volume—remember “green is good” and the circle will fill as you approach your target volume
3. Target ventilation rate defaults to provide feedback for a 30:2 compressions to ventilation scenario; **ALWAYS adjust** the monitor to 10 breaths/minute whenever advanced airway established AND providing continuous CPR
4. Place the monitor so the chest compressor and ventilation provider are in full view of monitor and able to adjust ventilation volume and rate accordingly
5. Remove ITD if ROSC is achieved; Re-apply if ROSC is lost
6. [Practical skill sheet](#)



Ben Weston, MD, MPH
Medical Director