



Burette / Buretrol Tubing


Indicated for blood administration in pediatric patients <40kg

Paramedic

PROCEDURE

- Ensure *all* clamps are closed
- Connect Buretrol spike into blood product
- Insert blood tubing spike into distal port of Buretrol
- Attach distal end of blood tubing to QinFlow Warrior Compact Disposable Unit (CDU)
- Prime blood tubing with saline
- Open proximal roller clamp to fill chamber with appropriate dosed blood volume
- Close proximal roller clamp; Open proximal vent
- Open Buretrol distal slide clamp to administer all blood products in fill chamber
- Repeat process until blood product dosing is complete

REFERENCE GRAPHICS

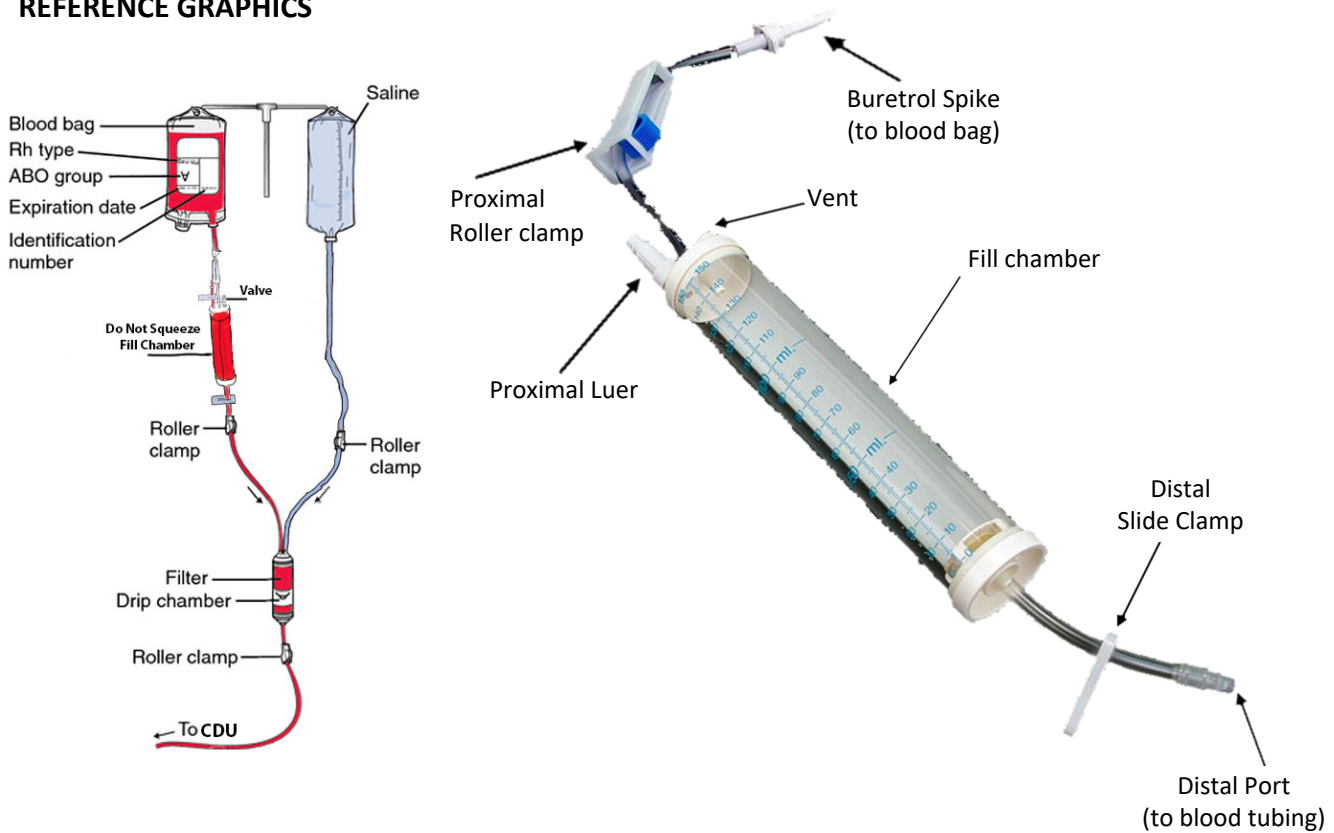

KEY POINTS

- Buretrols are indicated for blood administration in **pediatric patients <40kg**
- Not for use with other IV drips
- **Do not use pressure bag with Buretrol**
- **Do not squeeze the fill chamber**

Medical Director: Ben Weston, MD, MPH
November 2024